

SmartStand

**RAININ**  
Pipetting 360°

**The Intelligent Pipette Stand**

Check Calibration Status

View Service Data

Manage Compliance



**Feature Rich, Easy to Use**  
Pipette Management Simplified

**METTLER TOLEDO**

# Pipette Management Simplified


## SmartStand Fits Your Workflow

The Rainin SmartStand is an intelligent, easy-to-use pipette asset management system for busy labs and lab managers who want to:

- Reduce out-of-service/specification-related costs
- Proactively manage compliance
- Save time on recordkeeping and reporting

By reading the RFID chip that's in every Rainin XLS and XLS+ brand pipette, SmartStand instantly displays the compliance status for up to four pipettes – manual or electronic. Pick up a pipette and SmartStand's color LCD screen switches to a detailed view of that pipette's service and calibration status.

SmartStand works independently or integrates with EasyDirect™ Pipette Asset Management software to create a laboratory-wide pipette asset management and tracking system. EasyDirect simplifies lab asset management by enabling users to record pipette usage, flag pipettes for service, transfer data between other lab management systems (LIMS) and by recording unique data (such as owner, quick check methods and inventory number) onto each pipette's RFID chip.

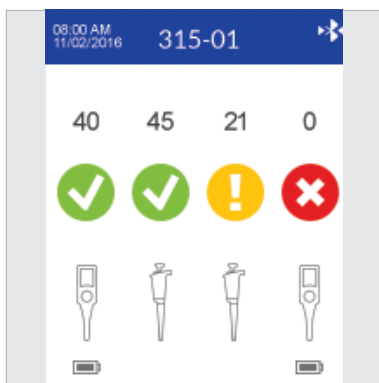

### Reduce Costs

Don't be caught off guard by pipettes that are out-of-specification or beyond their service dates. SmartStand brings clarity and transparency right to the benchtop.



### Manage Compliance

Transform tracking with an automated, highly-customizable system for proactively managing pipette service, calibration and periodic performance checking.



### Save Time

SmartStand simplifies and streamlines pipette management. It displays the service data that's stored on every Rainin XLS pipette and integrates with EasyDirect for a more comprehensive asset tracking and reporting system.



### SmartStand and Hang-Up in One!

SmartStand turns into a "Smart Hang-Up" with an accessory fastening kit.



# SmartStand and EasyDirect™

## Powerful Asset Management Tools

### SmartStand

- Flexible, scalable
- For manual and electronic XLS and XLS+ pipettes
- User defined service/calibration intervals
- QuickCheck reminders for periodic performance verification
- Intuitive graphical interface
- 110/220V, 50-60 Hz
- Connects via Bluetooth\* or USB
- Stands alone or use optional kits to mount the head to a shelf or wall

### EasyDirect Pipette Asset Management

- Connects up to multiple SmartStands to a PC
- Tracks and manages all brands/models of pipettes
- Reads/writes to RFID chips on Rainin XLS and XLS+ pipettes
- Export to Excel, csv or xml (append, overwrite or delete)
- Automatic export scheduler for autonomous operation
- Schedule service requests
- System requirements:  
Windows 7 (SP1) or higher, 5GB of free memory and 4 GB of RAM

Order No.	Description
30312897	Smart Charging Stand SCS-B
30312898	Smart Charging Stand SCS
30312899	Fastening Kit Magnetic SCS-FK-M
30313000	Fastening Kit Shelf SCS-FK-S
30313247	Fastening Kit Bracket SCS-FK-B
17012879	Elect Pipettes Pwr Supply E4-RCSWPS
30086494	ADP-BT-S Bluetooth RS232 Dongle

\* For SCS use Bluetooth dongle ADP-BT-S (Order No. 30086494)

[www.mt.com/rainin](http://www.mt.com/rainin)

For more information

### Mettler-Toledo Rainin, LLC

7500 Edgewater Drive, Oakland, CA 94621  
Phone +1 510 564 1600  
Fax +1 510 564 1604

Subject to technical changes  
© 10/2016 Mettler-Toledo Rainin, LLC  
Printed on demand from online content. PB-269 Rev A

